




<b>Product</b>	Showering and changing table 2000, coated canvas, electrically height adjustable
<b>Item number</b>	R8582418000 VVS:777544580 EAN:5708590355832
<b>Colour</b>	White
<b>Product description</b>	<p>Showering and changing table 2000, powered, height adjustable 700 mm; from 300 mm -1000 mm, with wired hand control Foldable. Incl. safety rail and water collection tray with integrated outlet. Outer measurements 900 mm x 1800 mm. Excl. fixing kit. See Technical datasheet (pdf) for information on fixings suitable for various wall types.</p> 
<b>Adjustability</b>	Height adjustable
<b>Adjustment height</b>	700 mm
<b>Product capacity</b>	<p>Maximum load: 175 kg. When the table is loaded with 175 kg, each mounting bolt supports 23 kg <i>This product has been tested statically with a load of 440 kg as the IEC 60601-standard requires testing with safety factor 2.5.</i></p>
<b>Materials</b>	Powder-laquered stainless steel, canvas (PVC-coated for extra durability) and polystyrene.
<b>Ergonomy</b>	<p>Pressalit Care offers a range of high quality, fixed or power adjusted showering and changing tables designed to provide a stable platform suitable for changing adults and children in a variety of settings: schools, day care centres, institutions, hospitals and private homes.</p> <p>The 2000 series nursing bench has a low start position, to make transfers to the bench easier, and with washable coated canvas for easy cleaning. The flexibility of the nursing bench provides for optimal working conditions for the carer as well as a high level of comfort and safety for the user.</p>
<b>Features</b>	<p>A water collection tray with integrated water outlet prevents water from splashing on carers or the floor. The flexible hose attached to the centre mounted outlet of the water collection tray, can be connected to either an outlet in the floor or on the wall. The integrated gas cylinders enable the bench to be folded out of the way when not in use to maximise floor space in front of the unit.</p> <p>Electric with wired hand control – the red buttons clearly indicate “Up” and “Down”. The wired hand control comes with a holder for fitting to the wall.</p> <p>The stretcher of the bench is very easy to clean, and the stretcher cover can easily be removed for washing or easy replacement.</p>
<b>Mains connection (V)</b>	230



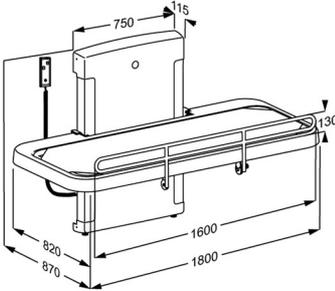






<b>Product</b>	Showering and changing table 2000, coated canvas, electrically height adjustable
<b>Power consumption (maximum load)</b>	Power consumption by maximum load: 70 W
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. <b>NOTICE</b> If the product is installed in an atmosphere heavy in salt or chloride, the warranty of the nursing bench will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product. The damaging effects can, however, be reduced if the product is showered with tap water on a daily basis.
<b>Cleaning</b>	<p>Showering and changing tables 1000 &amp; 2000. The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Clinimax Difficil-S Erisan OXY+</p> <p><b>NON-DISINFECTANTS:</b> Suma Star free D1 Taski 100 Taski 200</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>
<b>Declaration of conformity</b>	This product complies with the requirements of Council Directive 93/42/EEC concerning medical devices, as last amended by Council Directive 2007/47/EC of 5 September 2007. It also complies with all relevant requirements in the Council Directive 2006/42/EC (of 17 May 2006) on machinery and the Council Directive 2011/65/EC on RoHS.










<b>Product</b>	Showering and changing table 2000, coated canvas, electrically height adjustable
----------------	--

## Installation in different wall types

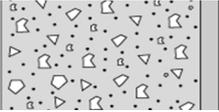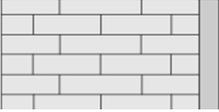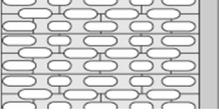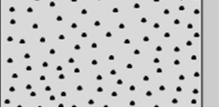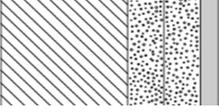	<b>Concrete &gt;B25</b> V8608 - Fixing kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.
	<b>Solid Brick &gt;= Mz12</b> V8608 - Fixing kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.
	<b>Solid Sand Lime Brick &gt;= Hlz12</b> V8608 - Fixing kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.
	<b>Aerated Concrete P4</b> V8658 - Fixing kit, 6 coach screws Ø10 x 70 - plugs GB14. To be ordered separately.
	<b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</b> V8608 - Fixing kit, 6 coach screws Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.



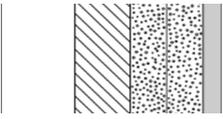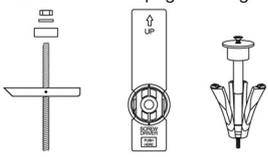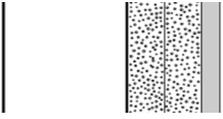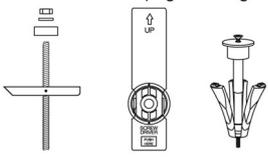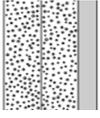

<b>Product</b>	Showering and changing table 2000, coated canvas, electrically height adjustable	
	<b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</b> Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.	fischer KD 8    snaptight    wingit 
	<b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm)</b> Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.	fischer KD 8    snaptight    wingit 
	<b>Drywall (2 x 12,5 mm) water-resistant</b> Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.	

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

